


PRECISION TWIN BLADE CUTTERS “DBC”

(Manually operated lever)

Designed to Paper, Foils, Films... samples, quickly and with just one operation of a certain width and 260mm length.
With a fixed blade and two movable cutter blades, made of special steel.

According to Standards: ISO 1924 – SCAN P38 – TAPPI...

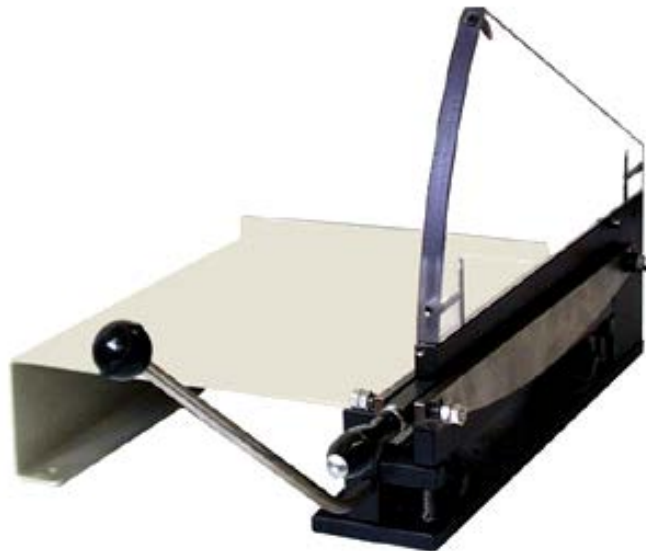

Before cutting, press the lever which meets the fixed cutter blade at the guide end. This clamps and tightens the sample and therefore ensures that high accuracy is obtained.

Equipment is **provided with methacrylate security cover protection**. This avoids the operator inserting his hands between the blades and the cutting base during cutting.

- ✓ DBC-12,7 model of **12,7 mm** cutting wide
- ✓ DBC-15 model of **15 mm** cutting wide
- ✓ DBC-20 model of **20 mm** cutting wide
- ✓ DBC-25 model of **25 mm** cutting wide
- ✓ DBC-25,4 model of **25,4 mm** cutting wide
- ✓ DBC-50 model of **50 mm** cutting wide

WEIGHTS AND DIMENSIONS:

Dimensions: 455 x 285 x 270 mm (W X D X H)
Box for transport: 600 x 450 x 450 mm (W X D X H)
Net/Gross Weight: 10 Kg /21 Kg

DELIVERY CONTENT:

> Twin Blade Cutter DBC... (model ordered)



* TECHLAB SYSTEMS reserves the right to do any technical modification without advance notice

Doc. : DBC_CUTTERS-1-CAT-E-R2