

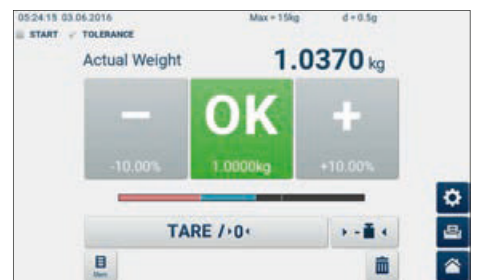


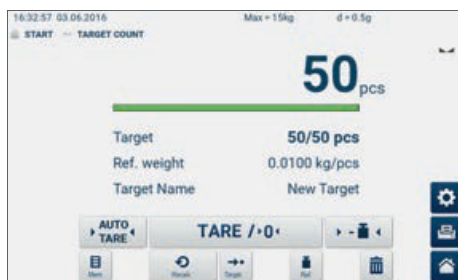

KERN Easy-Touch bench scale – the intuitive way to weigh

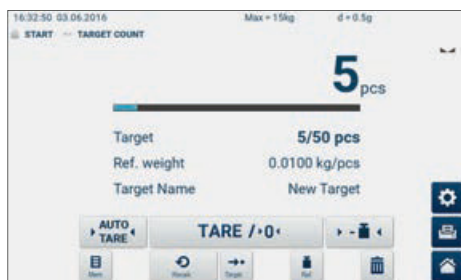

Convenient piece-counting function

Standard use: Direct input of the reference quantity or the reference weight.
 Professional use: Recall of items to be counted from the database with all relevant additional data such as piece weight, name, reference quantity, tare container weight, tolerances. Which means a super-fast sequence when counting: Select item – Load – Count – done!

Fill-to-target function: Programmable target quantity or target weight. A signal will be displayed when the target value is reached. A colour bar graph shows whether the quantity removed is below the limit (blue), is within the specified tolerances (green) or is above the limit (red). This means that the same target value can always be counted and picked, efficiently and accurately.

Convenient weighing with tolerance range (Checkweighing)

Standard use: Direct input of the tolerances in grammes or percent.
 Professional use: Recall of items from the database with all relevant additional data such as piece weight, name, reference quantity, tare container weight, tolerances. Which means a super-fast sequence when portioning, dosing or sorting: Select item – Load – Check – done!

Bench scale KERN GAT



Features

- Thanks to the **intuitive touchscreen concept** and multi-lingual operation (DE, GB, IT, ES, PT, NL, FI, PL, RUS, SE, CZ) the balance can be used easily by inexperienced users, straightaway. For the professional user, the balance provides convenient functions which permit a high level of individualisation and which thereby make operation significantly easier and quicker.
- Through a **large memory**, e.g. for item master data, weighing data etc. the balance is the ideal solution when working with a large range of goods or where there are high requirements for record keeping and documentation
- **1** Thanks to the high level of **connectivity** it is easy to connect barcode readers for item recognition, USB sticks for data storage or data transfer to PCs and networks using RS-232 or USB. This means that this range can be used for many different functions in Industry 4.0 applications
- **Convenient recipe weighing**
Standard use: Direct input of the recipe ingredients in grams or percent
Professional use: Recall of items from the database with all relevant additional data such as target value, tolerances, name, tare container weight. That leads to a super-fast sequence when recipe-weighing: Select recipe – Weigh – done!

- A graphic signal helps when dosing individual recipe ingredients. An information column on the lefthand edge of the screen allows you to see at a glance which ingredients have already been weighed (proportionately).
- **Take-Off function:** allows you to remove the same quantity every time, e.g. when portioning in cafeterias or when picking small parts. A colour bar graph shows whether the quantity removed is below the limit or within the specified tolerances (blue) or is above the limit (red)
- **Classifying:** Similar objects are automatically classified according to their weight into predetermined classes.
- **AUTO-PRINT function:** Automatic printing of the weighing result when the weighing value remains unchanged. The printing function as well as the content of the print-out can be adjusted individually depending on the operating mode
- **Remote control of the balance** using external switches or computers through KCP (KERN Communication Protocol), compatible with the large SICS command set (Standard Interface Command Set, Level 0/1)
- **2** The battery compartment of the device can be opened easily, which means it is possible to change the rechargeable batteries without having to use any tool or break any verification seal labels etc.

Technical data

- Large backlit LCD touch display, digit height 12 mm, screen diagonal 7" (155×85 mm)
- Dimensions weighing surface, stainless steel, W×D 300×225 mm
- Overall dimensions W×D×H 315×350×120 mm


Accessories

- **Rechargeable battery pack internal**, operating time up to 10 h with backlight, charging time approx. 12 h, KERN GAB-A04
- **Y-cable** for parallel connection of two terminal devices to the RS-232 interface on the scale, e.g. signal lamp or barcode reader and printer, KERN CFS-A04
- **Software** for rapid and easy processing of the database at the PC (exporting, processing, importing), KERN SET-1.0
- Further details, plenty of further accessories and suitable printers see *Accessories*

STANDARD





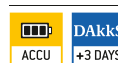






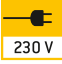





























OPTION

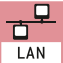

Model	Weighing range [Max] kg	Readout [d] g	Reproducibility g	Linearity g	Smallest part weight [counting] g/piece	Option	
						DAkkS Calibr. Certificate	DKD KERN
KERN							
GAT 6K-4	6	0,2	0,4	± 0,6	0,2	963-128	
GAT 10K-4	15	0,5	0,5	± 1,5	0,5	963-128	
GAT 30K-3	30	1	2	± 3	1	963-128	

KERN Pictograms:

 Internal adjusting: Quick setting up of the balance's accuracy with internal adjusting weight (motordriven).	 Piece counting: Reference quantities selectable. Display can be switched from piece to weight.	 Rechargeable battery pack: Rechargeable set.
 Adjusting program CAL: For quick setting up of the balance's accuracy. External adjusting weight required.	 Recipe level A: Separate memory for the weight of the tare container and the recipe ingredients (net total).	 Universal mains adapter: with universal input and optional input socket adapters for A) EU, GB B) EU, GB, CH, USA C) EU, GB, CH, USA, AUS
 Memory: Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.	 Recipe level B: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display.	 Mains adapter: 230V/50Hz in standard version for EU. On request GB, USA or AUS version available.
 Alibi memory: Electronic archiving of weighing results, complying with the 2014/31/EU standard.	 Recipe level C: Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display, multiplier function, adjustment of recipe when dosages are exceeded or barcode recognition.	 Power supply: Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request.
 Data interface RS-232: To connect the balance to a printer, PC or network.	 Totalising level A: The weights of similar items can be added together and the total can be printed out.	 Weighing principle: Tuning fork A resonating body is electromagnetically excited, causing it to oscillate.
 RS-485 data interface: To connect the balance to a printer, PC or other peripherals. High tolerance against electromagnetic disturbance.	 Percentage determination: Determining the deviation in % from the target value (100 %).	 Weighing principle: Electromagnetic force compensation Coil inside a permanent magnet. For the most accurate weighings.
 USB data interface: To connect the balance to a printer, PC or other peripherals.	 Weighing units: Can be switched to e.g. non-metric units at the touch of a key. See balance model. Please refer to KERN's website for more details.	 Weighing principle: Single cell technology Advanced version of the force compensation principle with the highest level of precision.
 Bluetooth* data interface: To transfer data from the balance to a printer, PC or other peripherals.	 Weighing with tolerance range: Upper and lower limiting values can be programmed individually for e.g. dosing, sorting and portioning.	 Verification possible: The time required for verification is specified in the pictogram.
 WLAN data interface: To transfer data from the balance to a printer, PC or other peripherals.	 Hold function: (Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value.	 DAkKS calibration possible (DKD): The time required for DAkKS calibration is shown in days in the pictogram.
 Control outputs (optocoupler, digital I/O): To connect relays, signal lamps, valves, etc.	 Protection against dust and water splashes IPxx: The type of protection is shown in the pictogram.	 Package shipment: The time required for internal shipping preparations is shown in days in the pictogram.
 Interface for second balance: For direct connection of a second balance.	 ATEX explosion protection: Suitable for use in hazardous industrial environments, in which there is explosion danger. The ATEX marking is specified for each device.	 Pallet shipment: The time required for internal shipping preparations is shown in days in the pictogram.
 Network interface: For connecting the scale to an Ethernet network. With KERN products you can use a universal RS-232/LAN converter.	 Stainless steel: The balance is protected against corrosion.	 Warranty: The warranty period is shown in the pictogram.
 Wireless data transfer: between the weighing unit and the evaluation unit using an integrated radio module.	 Suspended weighing: Load support with hook on the underside of the balance.	
 GLP/ISO log: The balance displays the weight, date and time, regardless of a printer connection.	 Battery operation: Ready for battery operation. The battery type is specified for each device.	
 GLP/ISO log: With weight, date and time. Only with KERN printers.		

KERN – Precision is our business

To ensure the high precision of your balance KERN offers you the the appropriate test weight in the international OIML error limit classes E1-M3 from 1 mg - 2500 kg. In combination with a DAkKS calibration certificate the best pre-requisite for proper balance calibration.

The KERN DAkKS calibration laboratory today is one of the most modern and best-equipped DAkKS calibration laboratories for balances, test weights and force-measurement in Europe.

Thanks to the high level of automation, we can carry out DAkKS calibration of balances, test weights and force-measuring devices 24 hours a day, 7 days a week.

Range of services:

- DAkKS calibration of balances with a maximum load of up to 50 t
- DAkKS calibration of weights in the range of 1 mg – 2500 kg
- Volume determination and measuring of magnetic susceptibility (magnetic characteristics) for test weights
- Database supported management of checking equipment and reminder service
- Calibration of force-measuring devices
- DAkKS calibration certificates in the following languages DE, GB, FR, IT, ES, NL, PL
- Conformity evaluation and reverification of balances and test weights

Your KERN specialist dealer: